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Canine/Feline Cognitive Dysfunction Rapid Test Kit

Adjunct to assessment of canine (CCDS) and feline (FCDS) cognitive dysfunction syndromes

Cat.no: 0082811

Colloidal gold double-antibody sandwich assay · Qualitative

Canine/feline serum or EDTA plasma · One C/T window · One S well · 1, 5, 10, 20 or 25 tests/box

Proprietary Technology

This IFU refers to the component recognized by the assay as the detection marker. Its molecular identity is proprietary.

Intended Use

For qualitative detection of a marker in canine or feline serum or EDTA plasma, as an adjunct to veterinary assessment of dogs and cats with cognitive or behavioral changes. Interpret alongside the history, behavioral assessment, physical examination and further testing as indicated.

Clinical Background

Cognitive dysfunction may affect orientation, the sleep–wake cycle, elimination habits, social interactions and activity. Consider other causes of similar signs, including pain, sensory decline, systemic disease and neurologic disorders.

Test Principle

The marker binds to a colloidal gold-labeled antibody and migrates along the membrane. A capture antibody at T binds the complex to form a visible colored line. The C line monitors flow through the membrane.

Kit Contents and Materials Required

Test cassettes, matched buffer and IFU are supplied. Check the pack for the exact configuration. Suitable droppers, collection and separation supplies, a timer and personal protective equipment are also required.

Storage and Expiry

Follow the storage conditions on the outer pack. Keep the pouch sealed and dry until use. Check the lot, expiry date and package integrity before testing.

Specimen Collection and Handling

Trained personnel should collect venous blood. For serum, allow the blood to clot and separate the serum. For plasma, collect in an EDTA tube, mix gently and separate the plasma. Transfer the clear fraction without disturbing cells or fibrin.

Record species, age, sampling time and relevant history. Use fresh specimens where possible. Recollect specimens that are clotted, contaminated, markedly hemolyzed or severely lipemic. For delayed testing, follow storage and transport conditions validated for this product.

Behavioral History to Record

Dogs: disorientation or getting stuck in familiar places, staring, altered sleep–wake cycles, nighttime pacing or vocalizing, house-soiling, altered interaction or anxiety, and reduced response to familiar commands.

Cats: unexplained nighttime vocalizing or pacing, increased daytime sleep, disorientation or staring in familiar places, soiling outside the litter box, reduced interaction or play, or anxiety.

Clinical Comparison

The case groups included animals with multiple owner-reported behavioral changes and a veterinarian's clinical assessment consistent with cognitive dysfunction. Healthy controls were 2–5 years old.

Dogs (CCDS)

Rapid test	Case group	Healthy controls	Total
Positive	27	0	27
Negative	1	50	51
Total	28	50	78

Positive agreement 27/28 · 96.43% (82.29–99.37%)
Negative agreement 50/50 · 100.00% (92.87–100.00%)
Overall agreement 77/78 · 98.72% (93.09–99.77%)

Cats (FCDS)

Rapid test	Case group	Healthy controls	Total
Positive	30	0	30
Negative	1	40	41
Total	31	40	71

Positive agreement 30/31 · 96.77% (83.81–99.43%)
Negative agreement 40/40 · 100.00% (91.24–100.00%)
Overall agreement 70/71 · 98.59% (92.44–99.75%)

Parentheses show two-sided Wilson 95% confidence intervals. Cases were classified by clinical assessment. A standardized reference protocol was not provided. The healthy controls were younger. These estimates describe the study groups and do not establish diagnostic performance in age-matched animals, animals with similar signs, or the general clinic population.

Zhejiang Anji Saianfu Biotech Co., Ltd.
2nd Floor, No. 3 Factory, No. 489 Wenyun Road, Tangpu Industrial Park, Dipu Subdistrict, Anji County, Huzhou City, Zhejiang Province, 313303, P.R. China. | www.safheal.com

Clinical Background References · Olby et al., JAVMA 264(4):498–505 (2026) · Ray et al., J Feline Med Surg 23(7):613–638 (2021)

Canine/Feline Cognitive Dysfunction Rapid Test Kit

Cat.no: 0082811

Colloidal gold double-antibody sandwich assay · Qualitative

Canine/feline serum or EDTA plasma · One C/T window · One S well · 1, 5, 10, 20 or 25 tests/box

Test Procedure

1. Verify specimen identity, expiry and package integrity. Bring the specimen, cassette and buffer to room temperature. Open the pouch and lay the cassette flat.
2. Use a dropper suitable for this kit to dispense 1 drop of serum or EDTA plasma into the sole S well. Avoid air bubbles.
3. Hold the buffer bottle upright. Add 3 drops of matched buffer to the same S well and start timing immediately after the last drop.
4. Keep the cassette level. In good light at 10 minutes, check C first, then read T. Results at 15 minutes or later are invalid.

Procedure Overview


Serum and EDTA plasma follow the same application steps.

01 Collect and separate

Use separated serum or plasma

02 Add specimen

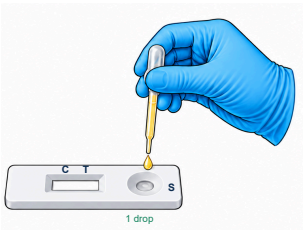1 drop specimen → S well

03 Add buffer

3 drops matched buffer → same S well

04 Time and read

10 min

≥15 min


✕

Read at 10 min · Invalid at ≥15 min

Result Interpretation

C = Control | T = Test

Positive · Marker detected

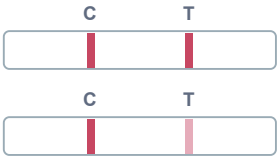Both C and T appear. Any clearly visible T line, including a faint line, is positive. With a compatible history and behavior, use the finding in further veterinary assessment for cognitive dysfunction.

Negative · Marker not detected

Only C appears. Persistent or progressive behavioral changes still warrant veterinary assessment.

Invalid

If C is absent, the result is invalid whether or not T appears. Investigate and repeat with a new cassette.

Illustrations are schematic. If background color obscures the result or the lines are unreadable, repeat with a new cassette.

Quality Control

C must appear for each test. It indicates migration through the membrane. Assess specimen quality separately. Run external controls as required by laboratory procedures.

Precautions

For trained personnel. Single-use cassette. Use the matched buffer. Handle and dispose of specimens and test materials as potentially infectious.

Limitations

This test is qualitative. Line intensity does not indicate marker concentration or disease severity. Results cannot independently confirm, exclude or stage cognitive dysfunction. If results and clinical findings differ, evaluate further and recollect and retest as appropriate.

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